

Compact Bench Scales

SJ-HS

Series



AND
A&D Australasia
"Excellence in Measurement"
www.andaustralasia.com.au

Features

- Large, removable stainless steel (SUS304) weighing pan, which is easy to wash and keep clean
- High-speed weighing in less than 1 second to increase your work efficiency
- Resolution up to 1/5000 (SJ-5001HS) for highly precise measurements
- Clearly visible LCD (with 25 mm character height) ensuring accurate reading of results
- Low-profile construction for storage in narrow shelves, which also helps when lifting heavy objects to the weighing pan
- Tare function to exclude the weight of the container from the weighing result
- Long battery life of approx. 2000 hours (optional AC adapter is also available)
 - Auto power off function so as not to waste battery power
 - Battery indicator that tells you the remaining battery level
- Four weighing units (g, lb, oz, and lb-oz) that can be switched with a single key
- Two side handles allowing easy carrying of the scale
- Zero tracking function to maintain the display at zero before measurement




The battery box makes it quick and easy to replace batteries.

Specifications (Trade Approved Models)

		SJ-1000HS	SJ-2000HS	SJ-5000HS	SJ-12KHS
Capacity x Resolution	g	1000g x 1g	2000g x 2g	5000g x 5g	10kg x 10g/12kg x 20g
Linearity		±1g	±2g	±5g	±10g
Repeatability (S.D.)		1g	2g	5g	10g
Span drift		±0.015% / °C TYP (5°C~35°C/ 41°F~95°F)			
Display		Liquid Crystal Display , 25mm / 0.98 inches			
Operating temperature		-10°C~40°C / 14°F~104°F, less than 85% RH (No condensation)			
Display update		Approx. 10 times per second			
Power		4 x R20P / LR20 / "D" size batteries or AC adapter			
Battery life		Approx. 800 hours with alkaline cells at 20°C / 68°F			
Pan size		230 (W) x 190 (D) mm / 9.05 (W) x 7.48 (D) inches			
Dimensions		265 (W) x 250 (D) x 103 (H) mm / 10.4 (W) x 9.8 (D) x 4.1 (H) inches			
Weight		Approx. 2.1kg / 4.6lb			
Calibration mass		1000g	2000g	5000g	10kg

Specifications subject to change for improvement without notice

External Dimensions (mm / inches)

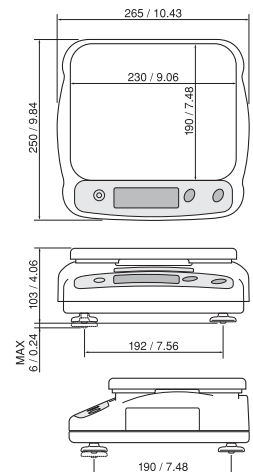

Option AC adapter



"Excellence in Measurement"

ACN 26 007 556 809

Head Office
32 Dew Street
THEBARTON South Australia 5031
Telephone (08) 8301 8100
Facsimile (08) 8352 7409

Victorian Office
39 Bakehouse Road
KENSINGTON Victoria 3031
Telephone (03) 9371 1555
Facsimile (03) 9372 1193

New South Wales Office
Unit 4, 14 Abbott Road
SEVEN HILLS New South Wales 2147
Telephone (02) 9674 5466
Facsimile (02) 9674 2544

sales@andweighing.com.au
www.andweighing.com.au

Your Dealer